



SEM

SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724

Mobile 1: +91 98791 03905

Email: sem@semtrainers.com

Website: www.semtrainers.com


IV Arm for 1005650

Item No. LF03214U
Weight 5.44kg
Make: Nasco Healthcare, USA

[Read More](#)

SKU:

Categories:ALS Adult

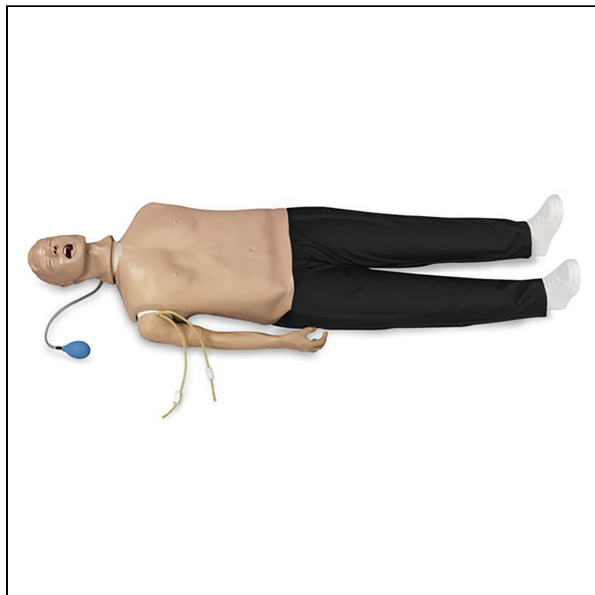

Product Description

The IV Arm is available in a version that connects to the right arm socket of the CPARLENE® manikin. Features realistic vinyl skin and latex veins so the skin actually rolls as you palpate the vein, and the characteristic “pop” can be felt as the needle penetrates the vein. Realistic flashback occurs with proper insertion. The arm pivots at the elbow for easy accessibility to antecubital fossa, along the forearm, and at the back of the hand. A bony landmark at the shoulder helps identify muscle tissue for intramuscular injections. Replacement skins and vein sealant are available. The kit includes the Injectable Training Arm with white skin, two fluid supply bags, 3 cc syringe, 12 cc syringe, needle, one pint blood, teaching guide, and carrying case.

Life/form® IV Arm for CPARLENE® Manikin 1005650 (Item No. LF03214U | Weight: 5.44 kg | Nasco

Healthcare, USA)

A complete intravenous training arm module designed to connect directly to the right arm socket of the CPARLENE® manikin (model 1005650). This arm delivers high-fidelity IV insertion, blood draw, and intramuscular injection practice within the context of a full-body manikin simulation — integrating vascular access training into complete patient care scenarios rather than isolating it to a standalone task trainer.

Simulation Fidelity:

The arm features realistic vinyl skin and latex veins so the skin actually rolls as you palpate the vein, and the characteristic "pop" can be felt as the needle penetrates the vein. Realistic flashback occurs with proper insertion. This tactile feedback is what separates competency-building from mere procedural repetition — students develop the palpation sensitivity and needle control required on real patients. [semtrainers](#)

Access Points & Anatomical Features:

- The arm pivots at the elbow for easy accessibility to the antecubital fossa, along the forearm, and at the back of the hand [semtrainers](#)
- A bony landmark at the shoulder helps identify muscle tissue for intramuscular injections [semtrainers](#)
- Multiple venous access sites support repeated practice across antecubital, forearm, and dorsal hand locations

Kit Contents:

Includes the Injectable Training Arm with white skin, two fluid supply bags, 3 cc syringe, 12 cc syringe, needle, one pint simulated blood, teaching guide, and carrying case. [semtrainers](#)

Consumables Available Separately:

Replacement skins and vein sealant are available — essential for maintaining simulation fidelity after high-frequency use. See the Injectable Training Arm Replacement Skin and Vein Kit for compatible consumables. [semtrainers](#)

Ideal for: Nursing colleges, paramedic training programmes, MBBS clinical skills labs, and hospital in-service training departments requiring integrated IV and IM injection practice within a full-body manikin setup.

For institutional or bulk orders, contact SEM Trainers.