



SEM

SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724
Mobile 1: +91 98791 03905

Email: sem@semtrainers.com
Website: www.semtrainers.com




Catheterization Simulator BASIC, female

Item No. 1020231 [P93B-F]
Weight 2.3kg
Dimensions 36 x 11 x 26.5 cm
Brand 3B Scientific

[Read More](#)

SKU:

Categories:Catheterization

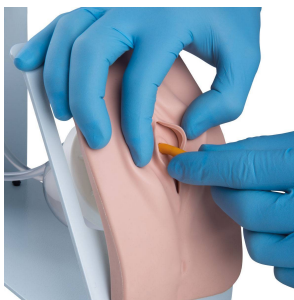

Product Description

With the female P93 BASIC catheterization trainer, female bladder catheterization procedures can be realistically demonstrated, practiced and assessed. The genital insert is placed into the holder and held in place with magnets. The material of the insert is soft and flexible: the labia are movable and can be spread. This way the student can practice all the necessary gestures for catheterization (e.g.

disinfection). The material gives a realistic feeling when inserting and removing the catheter into the bladder with realistic resistance and insertion depths. **Direct learning:** When the catheter has been correctly inserted, fluid runs out just like with a real patient. The catheter can also be checked through the transparent bladder at any time. **All the benefits of the female P93 BASIC catheterization simulator at a glance:**

- Female catheterization with realistic resistance
- Soft and movable labia
- Liquid outflow if catheterization is successfully carried out
- The situation can be checked through the transparent bladder
- Stands securely on the table thanks to the non-slip feet
- Magnetic connectors for quick set-up and dismantling
- Easy to clean, and the inserts can be completely removed

Delivery contents: • Stands with magnetic connections for genital inserts • Catheterization set with transparent bladder • Bottle and tube connectors **Note:** For catheterization, we recommend using catheters with only one or two drainage eyes opposite each other. Both silicone and PVC catheters can be used. We recommend the following catheters:

- Urethral narrowing “open” = 16 CH catheter or 14 CH balloon catheter
- Urethral narrowing “partial narrowing” = 14 CH catheter or 12 CH balloon catheter

Catheters with a curved tip should not be used! Tip: According to the guidelines for correct catheterization, even during the simulation, lubricating gel (not included) should be placed into the urethra to ensure that the catheter can be inserted smoothly. Alternatively, some lubricating gel can be placed on the tip of the catheter.

Clean, detailed, and clinically precise — good page, no data conflicts. Only thing worth flagging: preserve the catheter-sizing table exactly as given, since these are real clinical specification numbers, not marketing filler.

[Cert] Item No. 1020231 [P93B-F], weight 2.3 kg, dimensions 36×11×26.5 cm, Brand 3B Scientific. This is a bench-top (not full-body) female catheterization trainer using a magnetically-mounted, interchangeable genital insert with movable/spreadable labia, realistic insertion resistance and depth, and fluid outflow feedback on successful catheterization — visible confirmation through a transparent bladder.

[Cert] The source gives specific, actionable clinical guidance worth preserving verbatim rather than paraphrasing loosely: catheters with only one or two opposing drainage eyes are recommended (silicone or PVC), sized per urethral narrowing scenario — 16 CH catheter or 14 CH balloon catheter for “open,” 14 CH catheter or 12 CH balloon catheter for “partial narrowing” — and curved-tip catheters

should not be used. Lubricating gel is recommended but not included.

DESCRIPTION (copy-paste into WordPress):

Why Choose the Catheterization Simulator BASIC, Female, for Your Skills Lab in India?

The Catheterization Simulator BASIC, female (Item No. 1020231 [P93B-F], 2.3 kg, 36 × 11 × 26.5 cm, Brand: 3B Scientific) is a compact, table-top trainer for realistic demonstration, practice, and assessment of female bladder catheterization procedures. The genital insert mounts on the stand via magnetic connectors for quick setup and teardown, with soft, flexible, movable labia that let students practice the full procedural sequence — including disinfection — before catheter insertion.

The insert material replicates realistic insertion resistance and depth, and successful catheterization produces genuine fluid outflow, just as it would with a real patient — giving students immediate, unambiguous feedback on technique. Catheter placement can also be visually confirmed at any time through the transparent bladder.

Features:

- Female catheterization with realistic resistance
- Soft, movable labia
- Fluid outflow confirms successful catheterization
- Transparent bladder for visual confirmation
- Non-slip feet for secure table placement
- Magnetic connectors for quick setup/teardown
- Fully removable inserts for easy cleaning

Delivery contents: Stand with magnetic connections for genital inserts, catheterization set with transparent bladder, bottle and tube connectors.

Recommended catheters (not included):

- Urethral narrowing "open": 16 CH catheter or 14 CH balloon catheter
- Urethral narrowing "partial narrowing": 14 CH catheter or 12 CH balloon catheter
- Use catheters with only one or two opposing drainage eyes. Both silicone and PVC catheters are suitable. **Do not use catheters with a curved tip.**
- Lubricating gel is recommended (not included) for smooth insertion, applied to the urethra or the catheter tip.

Why Buy from SEM Trainers? SEM Trainers & Systems is an authorized distributor of 3B Scientific

products in India with over 30 years of experience supplying medical simulators and clinical skills training equipment to 1,500+ institutions across 28 states. We provide pan-India delivery and dedicated after-sales support.

Q: Which catheters should I use with this simulator? A: 16 CH or 14 CH balloon for open urethral narrowing scenarios, 14 CH or 12 CH balloon for partial narrowing. Avoid curved-tip catheters.

Q: Is lubricating gel included? A: No — lubricating gel is recommended for realistic practice but must be sourced separately.

Q: How do I order from SEM Trainers? A: Contact us at sem@semtrainers.com or call +91-88495 63724 for pricing and availability.