



SEM

SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724
Mobile 1: +91 98791 03905

Email: sem@semtrainers.com
Website: www.semtrainers.com





CVC Insertion Simulator II

Item No. M93UB
Weight 9 kg
Brand KYOTO KAGAKU, JAPAN

[Read More](#)

SKU:

Categories:Ultrasound Phantom

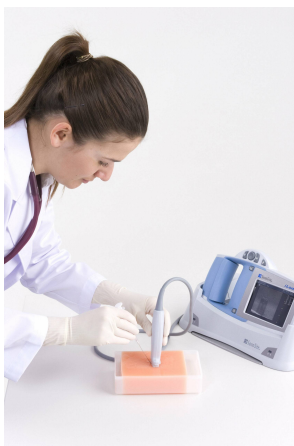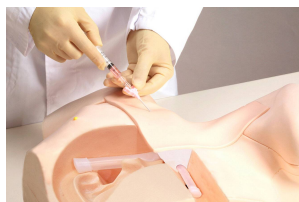

Product Description

Why Choose the CVC Insertion Simulator II for Your Vascular Access Training Program in India?

The CVC Insertion Simulator II (Item No. M93UB, 2.3kg product weight, Brand: Kyoto Kagaku, Japan) lets trainees conduct a full variety of central venous catheter insertion and management simulations using three interchangeable pads — covering both ultrasound-guided and landmark-guided technique on one platform.

The landmark puncture pad features anatomically correct vein bifurcations that simulate real mechanical complications — pneumothorax, catheter mislodging, and arterial puncture — a genuinely important safety-training feature, letting trainees learn to recognize and avoid these serious complications before ever encountering them with a real patient. The transparent anatomical block supports deeper anatomical understanding and guide wire manipulation practice.

Training applications: Ultrasound-guided CVC, landmark-guided CVC, ultrasound-guided venous access, prevention of mechanical complications.

Set includes: Male upper torso manikin, landmark puncture pad, ultrasound puncture pad, transparent anatomical block, introductory ultrasound training block, 10 red coloring cotton swabs, 10 blue coloring cotton swabs, air bulb, 2 plastic jars, irrigation bottle, syringe, sample syringe, sample needle.

Please note: an ultrasound scanner and needles are **not included**.

Specifications: Size W40×D20×H34cm (W15.8×D7.9×H13.4in). Product weight 2.3kg (5lb). Materials: soft resin, latex-free.

Replacement parts available: Landmark puncture pads (11347-119), ultrasound puncture pads (11347-170), introductory ultrasound training blocks (11347-210), irrigation bottle (11347-320), red coloring swabs (11388-400), blue coloring swabs (11402-010), pulsation unit for M93UB (11347-260).

Optional parts: Pulsation unit for CVC placement pad (11347-250), CVC placement pads (11347-240).

Why Buy from SEM Trainers? SEM Trainers & Systems is an authorized distributor of Kyoto Kagaku products in India with over 30 years of experience supplying medical simulators and vascular access training equipment to 1,500+ institutions across 28 states. We provide pan-India delivery and dedicated after-sales support.

Frequently Asked Questions

Q: What complications can be simulated? A: Pneumothorax, catheter mislodging, and arterial puncture — helping trainees learn to recognize and avoid these before working with real patients.

Q: Is an ultrasound scanner included? A: No — the scanner and needles are not included.

Q: How do I order from SEM Trainers? A: Contact us at sem@semtrainers.com or call +91-88495 63724 for pricing and availability.

Conduct a variety of CVC insertion and management simulations with the assistance of 3 interchangeable pads

Features

Training skills / Applications

Set includes

Size (approx.)

Packing size (approx.)

Weight (approx.)

Packing weight (approx.)

Materials

Replacement parts

Optional parts

Related products

Production ? Development Supervision

Notes

Update

1. Realistic ultrasound-guided CVC 2.Landmark puncture pad with anatomically correct vein bifurcations simulates mechanical complications including pneumothorax, mislodging and artery puncture 3.Transparent anatomical block for anatomical understanding and guide wire manipulation

Ultrasound-guided CVC Landmark guided CVC Ultrasound-guided venous access Prevention of mechanical complications

1 male upper torso manikin / 1 landmark puncture pad / 1 ultrasound puncture pad / 1 transparent anatomical block / 1 introductory ultrasound training block / 10 cotton swab for coloring (red) / 10 cotton swab for coloring (blue) / 1 air bulb / 2 plastic jars / 1 irrigation bottle / 1 syringe / 1 sample syringe / 1 sample needle

W40 x D20 x H 34cm / W15.8 x D7.9 x H 13.4in

W52x D45x H38cm/W20x D18x H15in

2.3kg / 5lb

9kg / 19.8lb

Soft resin / Latex free

11347-119 2 landmark puncture pads / 11347-170 2 ultrasound puncture pads / 11347-210 2 introductory ultrasound training blocks / 11347-320 irrigation bottle/11388-400 cotton swab for coloring (red, set of 10)/ 11402-010 cotton swab for coloring(blue, set of 10)/11347-260 1 pulsation unit (for M93UB)

11347-250 pulsation unit for CVC placement pad / 11347-240 2 CVC placement pads

[MW18 Ultrasound-Guided PICC Training Simulator M93C CVC Insertion Simulator III](#)

Masahiro Tanabe, M.D., Ph.D., Professor Director General Medical Education Center Chiba University School of Medicine / Kinya Sando, M.D., Ph.D., Professor Director Department of Human Dietetics Graduate School of Human Science Osaka Shoin Women's University / Masanori Hoki, M.D., Ph.D., Professor Department of Human Dietetics Graduate School of Human Science Osaka Shoin Women's University / Joho Tokumine, M.D.Ph.D., Department of Anesthesia, Kawatetsu Chiba Hospital

Ultrasound scanner and needles are not included.

October 1, 2020