




## SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

**Mobile :** +91 88495 63724  
**Mobile 1:** +91 98791 03905

**Email:** sem@semtrainers.com  
**Website:** www.semtrainers.com

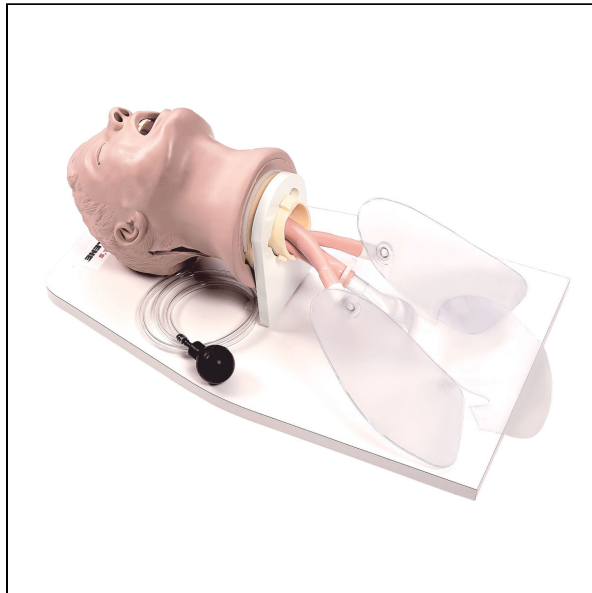

## Life/form® “Airway Larry” Adult Airway Management Trainer with Stand

Item No. LF03699U  
Weight 11.439kg  
Dimensions 63 x 30.5 x 40.6 cm  
Make: Nasco Healthcare, USA

[Read More](#)

**SKU:**

**Categories:**Airway Management Adult

## Product Description

The Life/form® Airway Management Trainer simulates non-anesthetized patients for practicing intubation, ventilation, and suction techniques. Realistic anatomy and landmarks including teeth, tongue, oral and nasal pharynx, larynx, epiglottis, arytenoids, false cords, true vocal cords, trachea, lungs, esophagus, and stomach. The trainer allows you to practice oral, digital, and nasal intubation, as well as E.T., E.O.A., P.T.L., L.M.A., E.G.T.A., and Combitube® insertion. Suction techniques and proper cuff inflation may also be practiced and evaluated. Features durable, rugged, one-piece construction and bifurcated lungs. “Airway Larry” represents an adult with more room to maneuver than the Adult Trainer (Number available upon request). With its slightly anterior position, “Airway Larry” is a great trainer for introductory as well as advanced students. It is ideal for testing purposes and Combitube® practice. Applying pressure to the cricoid cartilage changes the position of the trachea and closes the esophagus, allowing realistic practice of the Sellick maneuver. The jaw thrust technique can

also be practiced. The trainer can also be attached to Nasco's CPARLENE®, CRiSis™ Manikin. Comes with training stand to mount head, pump spray lubricant, and hand carrying hard case.

	Head/torso	Ideal for	Digital intubation	Laryngospasm & Tonge edema	Stomach rise	Sucti-oning	Jaw thrust	Manual carotid pulse	E.T., E.O.A., P.T.L, L.M.A., Combitube®, and King System insertion	Additional features
<b>Adult Airway Management Trainer Head</b> <i>Item No: <a href="#">1012740</a></i>	male head	introductory classes	Yes	-	-	Yes	-	Yes	Yes	can be attached to any of the adult CRiSis® Manikins
<b>Advanced Trainer Head</b> <i>Item No: <a href="#">1020961</a></i>	male head mounted on a board	introductory and advanced classes	Yes	Yes	Yes	Yes	-	-	Yes	can be attached to any of the adult CRiSis® Manikins
<b>Adult Airway Management Trainer</b> <i>Item No: <a href="#">1005633</a></i>	male head mounted on a board	testing purposes and Combitube® practice	Yes	-	-	Yes	Yes	-	Yes	can be attached to any of the adult CRiSis® Manikins
<b>Economy Adult Airway Management Trainer</b> <i>Item No: <a href="#">1005780</a></i>	male head mounted on a board	rescue & advanced airway management training	-	-	-	Yes	-	-	only LMAs and Combitube®	durable construction and bifurcated lung
<b>Adult Airway Management Trainer</b> <i>Item No: <a href="#">1017871</a></i>	torso with a female head	difficult airway and CPR training	Yes	-	-	Yes	Yes	Yes	Yes	durable construction and bifurcated lung
<b>Advanced Trainer Head</b> <i>Item No: <a href="#">1019249</a></i>	male head mounted on a board	advanced airway management training	Yes	Yes	-	Yes	-	Yes	Yes	durable construction and bifurcated lung

<b>S.A.L.A.D. Simulator</b> <i>Item No:</i> <a href="#">1021583</a>	male head	suction assisted laryngoscopy and airway decontamination	Yes	-	-	Yes	-	Yes	Yes	rigid suction catheter distracts lower mandible and tongue to permit laryngoscope insertion
<b>Airway Management Trainer</b> <i>Item No:</i> <a href="#">1019928</a>	male torso	airway management and CPR training	Yes	-	-	Yes	-	Yes	Yes	durable construction and bifurcated lung
<b>Adult Deluxe Airway Management Trainer</b> <i>Item No:</i> <a href="#">1018851</a>	male head mounted on a board	cricothyrotomy and advanced airway management	-	Yes	Yes	-	-	-	-	lifelike skin that will not be damaged by aggressive intubation
<b>Adult Airway Management Trainer</b> <i>Item No:</i> <a href="#">1005774</a>	male torso mounted on a board	rescue airway and anatomy training	-	Yes	Yes	-	-	-	-	review of the pulmonary system and oral/nasal pharyngeal spaces