



SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724
Mobile 1: +91 98791 03905

Email: sem@semtrainers.com
Website: www.semtrainers.com




Difficult Airway Management Simulator -Training Model-

Item No. MW13
Weight 16kg/35.2lb
Make: Kyoto Kagaku, Japan

[Read More](#)

SKU:

Categories: Airway Management Adult

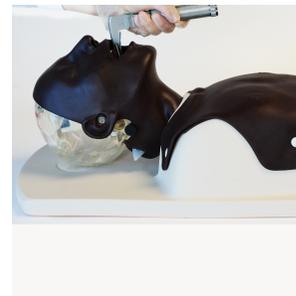

Product Description

From beginners to experts with a range of difficulty levels/Difficult Airway Management Simulator provides training for users

Features

1. 24 variations of patient scenario (including 1 normal case): 3 stages of mouth opening, 2 stages of neck flexibility, 2 tongue sizes and 2 positions of the vocal cords
2. Upper incisors are designed to break off when excessive force is applied.
3. For Nasal intubation: anatomically correct middle and inferior nasal concha to lean direction of the tube.

Training skills / Applications

Airway opening techniques (head tilt, jaw thrust) Bag-Valve-Mask ventilation Pre-intubation airway assessment Recreating the sniffing position Pressurization of external larynx to improve the laryngeal view Intraoral/Intranasal Intubation Use of oropharyngeal airway (OPA) Use of nasopharyngeal airway (NPA) Use of laryngeal airway mask Use of video laryngoscope Confirmation of successful ventilation by: -observation of thoracic and abdominal movement (lung expanding, stomach inflating) or auscultation of the chest Simulates incorrect procedures including intubation in the esophagus and unilateral intubation Practice securing tubes after intubation

Case / Pathology

Variation of DAM Setting Neck Flexibility (Life-like jaw movement) Mouth Opening/ - Normal /- Intermediate /- Difficult Tongue/ - Normal/- Swollen Laryngospasm/ - Normal/ - Laryngospasm

Set includes

1 manikin/ 3 upper incisors /2 tongues II/ 1 lubricant/ 1 syringe/ 1 carrying bag/ 1 instruction manual

Size (approx.)

W70 x D39 x H25 cm / W27.6 x D15.4 x H9.8 in

Packing size (approx.)

W74 x D46 x H47cm/W29 x D18 x H18.5in

Weight (approx.)

7.0kg / 15.4lb

Packing weight (approx.)

16kg/35.2lb

Materials

Soft resin/Latex free

Replacement parts

11392-010 10 upper incisors / 11392-090 1 face mask/11392-040 1 chest cover/11392-080 2 replacement tongues II / 11392-050 5 pair of lungs/11392-060 5 stomachs/11229-050 lubricant

Optional parts

11384-100 bronchofiberscopy training unit

Related products

[MW14 Difficult Airway Management Simulator -Demonstration Model-](#) [MW11 Difficult Airway Management Simulator -Assessment System-](#)

Recommended Devices

"Macintosh laryngoscope: blade size 3, 4 / Endotracheal tube: internal diameter 7.0mm, 7.5mm / Laryngeal mask Air-Q: size 3.5, LMA Supreme: size 4, I-gel: size 4"

Update

October 20, 2020