




SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724
Mobile 1: +91 98791 03905

Email: sem@semtrainers.com
Website: www.semtrainers.com

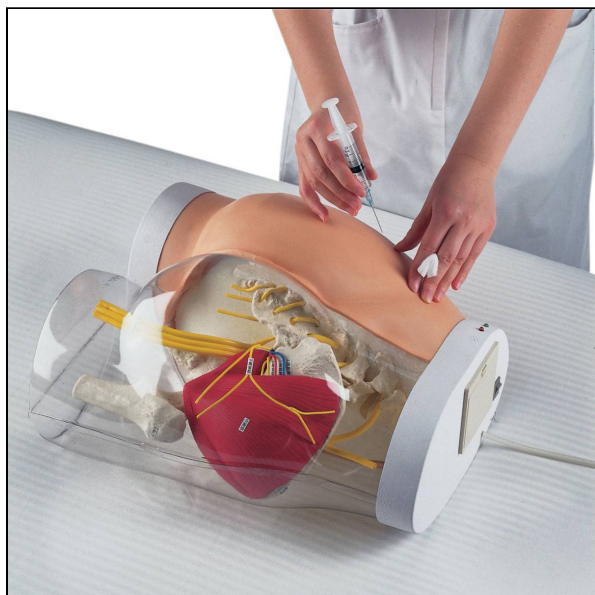


Two-in-One i.m. Injection Model of Buttock

Item No. 1005394 [W30503]
Weight 6.799 kg
Dimensions 38 x 35 x 22 cm
MPN: M152

[Read More](#)

SKU:

Categories: Injections and Punctures

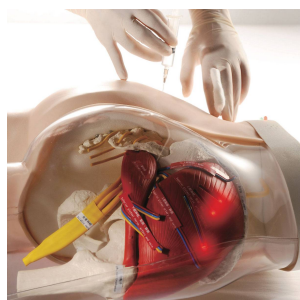

Product Description

This unique model of a human buttock has a special two-in-one function: On the right hand side, the anatomy of the bones, ilium crista, greater trochanter, m. gluteus medius, nerves and veins can be studied through the transparent outer structure. The student can transfer what he has learned about the position of nerves, veins etc. to the left hand side on which intra-muscular injections can be practiced. Anatomical landmarks can be palpated through the soft skin to identify the correct positions for injection. The injection of fluid (water) is possible. Correct injections are confirmed by a green light,

if the injection is placed in a wrong position or too deep, a red light will appear and a buzzer will sound.

	Intramuscular Injection Simulator - Upper Arm	Intramuscular Injection Simulator - Upper Leg	Buttock Injection Simulator	Intramuscular injection simulator	Intramuscular Injection Model of Buttocks	Intramuscular Injection Simulator	Intradermal Injection Simulator
Item No.	1009840	1000511	1000514	1010008	1005394	1005586	1005628
Body part	upper arm	upper leg	buttock	muscle extract with skin	pelvis with buttocks	lower torso	forearm
Palpable anatomical landmarks	✓	✓	✓	-	✓	✓	-
Injection type	intramuscular	intramuscular	intramuscular	intramuscular	intramuscular	intramuscular	intradermal
Injection sites	as indicated by a trainer	as indicated by a trainer	as indicated by a trainer	all the surface	as indicated by a trainer	as indicated by a trainer	eight
Proper injection	audiovisual feedback	audiovisual feedback	audiovisual feedback	-	green light will appear	-	skin welt will form
Additional information	training and test modes are available	training and test modes are available	training and test modes are available	will absorb fluid	special skin material does not show needle marks	simulated bony structure is embedded in the torso	liquid skin repair kit included (M-1013686)