



SEM

SEM TRAINERS & SYSTEMS

WHERE TIME MOVES AHEAD TO KEEP PACE WITH KNOWLEDGE

SEM- Scientific Educative Methods in Science, Engineering & Medicine

Mobile : +91 88495 63724

Mobile 1: +91 98791 03905

Email: sem@semtrainers.com

Website: www.semtrainers.com

Prestan Professional Adult Jaw Thrust CPR-AED Training Manikin Medium Skin

Item No. PP-JTM-100M-MS

Weight

Make: Prestan Products, USA

[Read More](#)

SKU:

Categories:AED Trainers



Product Description

The PRESTAN Professional Jaw Thrust Head gives instructors the option to teach with the jaw thrust maneuver or head-tilt/chin-lift allowing the student to grasp the manikin's lower jaw and lift with both hands displacing the jaw of the manikin forward and upward. The jaw thrust maneuver is the safest approach to opening the airway of a victim who has a possible head, neck or spinal injury because it can be accomplished by minimizing neck movement.

Jaw Thrust Heads may be purchased separately and easily fit onto the torso of the original PRESTAN

Professional Adult Manikin.

Fully compliant with current industry guidelines, including the most recent AHA Integrated Feedback Directive. **Prestan Professional Adult Jaw Thrust CPR-AED Training Manikin (installed with Rate Monitor and Depth Clicker for CPR Feedback and AA Batteries not included). Includes :**

- 10 Adult Face-Shield/Lung-Bags
- Instruction Sheet
- Carry Bag

Warranty Information & Guarantee


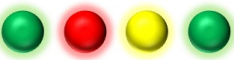






PRESTAN offers a THREE YEAR LIMITED WARRANTY against any defects in workmanship or materials under normal use with respect to its PRESTAN TRAINING MANIKINS.

Additional Details:


All PRESTAN Professional Manikins are available with a revolutionary CPR Rate Monitor that allows for instant feedback to both instructors and students. The visual feedback from the CPR Rate Monitor gives students a real life feel of delivering 100-120 compressions per minute. Lights continually signal rate progression, confirming compressions are within the 100- 120 cpm range recommended by the Industry Guidelines. The yellow light will blink if the 120 cpm upper limit is exceeded, warning students to slow down. Instructors can monitor several students quickly and easily.

In addition to the visual CPR Rate Monitor, students will also hear a clicker sound as the chest is pushed to the appropriate depth of 2.0 -2.4 inches allowing them to experience the true force needed to deliver real life chest compressions.

The CPR Rate Monitor may be purchased separately to upgrade any non-monitored PRESTAN Professional Manikin.

Sleep Mode/No Power	
1st Full Compression (Wake mode/LED Test Mode)	
CPM* < 60	
CPM* 60 – 79	
CPM* 80 – 99	
CPM* 100 - 120	
CPM* > 120	
Compressions Stopped	

*CPM = Compressions per minute

 = Flashing

